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USER'S GUIDE

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THANK YOU FOR YOUR PURCHASE
AND WELCOME TO NEOPULSE!

SPECIAL WARNINGS

DANGER

• Check the condition of the USB-C cable power cord before first use and then regularly during the life of the device. • Do not handle the USB-C cable with wet hands. • Do not immerse the USB-C cable or drop it into water. • If the toothbrush gets wet, dry it before plugging it into the power outlet to charge it.

WARNING


• Do not connect this appliance to a system with a different voltage different from that specified on the unit. • Do not use this product if the cord or plug is damaged. • Keep the USB-C cable away from hot surfaces. • NEOSONIC® toothbrushes from NEOPULSE® are comply with safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, consult your physician before using this device. • Use the appliance as explained in the instructions for use, this toothbrush has been designed for teeth and gums brushing only. Do not use it for any other purpose. • If you have had dental treatment in the last three months, consult your dentist, periodontist or pharmacist before using this toothbrush. • Give specific instructions to children (from 10 years old), as well as people with special needs regarding the proper use of this product and supervise them. • Do not drop the device. • This appliance does not contain any user-replaceable components. • If the product is opened or taken apart for any reason, the warranty will be invalidated. • Do not put the brush head, handle or charger in the dishwasher. • Replace the brush head at least every 3 months and earlier if you notice that the blue bristles of the brush whiten by 1 millimetre. • Use only NEOKIDS® or NEOSONIC® brush heads from NEOPULSE®.

KEEP THESE INSTRUCTIONS.

To stop your toothbrush, press the setting button (2) again after 5 seconds.

Step 3: Brushing

Apply a small amount of toothpaste. Position the bristles on the teeth at a slight angle (45 degrees) along the gum line and press the On/Off and Mode button (2).

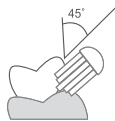For an even brushing divide your mouth into 4 distinct sections and spend 30 seconds brushing on each zone (4 sections X 30 seconds. = 2 minutes).

With Smart Brushing, every 30 seconds, a gentle vibration will tell you to change your brushing zone.

The NEOSONIC® toothbrush automatically stops after 2 minutes of brushing cycle, thus respecting the recommendations of dentists.

Brushing recommendations:

To clean the outer surfaces: brush top and bottom separately. Move the brush in a rotating motion, always from pink to white, to brush your teeth and gums without exerting pressure and especially do not rub.

To clean the inner surface of the teeth: make several vertical brushings on each tooth, always from the gum to the tooth. Continue for the duration of the timer (2 minutes), until the toothbrush stops. When using a sonic toothbrush for the first time, it is possible to at the beginning to feel some tickling. This feeling is common and disappears quickly after several uses.

Step 4: Our advice

Once the 2-minute brushing cycle is complete, you can go back and insist on the stain-prone areas, but also on the chewing surfaces of the teeth and in crannies where you have a tendency to get tartar.

Change your brush head at least every 3 months. When the blue bristles of the brush whiten over 1 millimetre, it is time to replace it.

NEOSONIC® toothbrushes are equipped with 3D brush heads with rounded bristles for efficiency in plaque removal while taking care of sensitivity of your teeth and gums. The colored bristles in the center serve as a visual guide for toothpaste dosing.


PRODUCT DESCRIPTION

1. Brush head with bristles wear indicator
2. 5 brushing modes & On / Off
3. 3 strength settings
4. Brushing mode indicator
5. Battery level indicator
6. USB-C charger plug
USB-C cap provided.
Optional use - the device remains fully functional without it.

Accessories:

- Hygienic protective cover
- 2 brush heads:
 - 1 soft (0.152 mm)
 - 1 ultra-soft (0.127 mm)
- 1 USB-C cable
- 1 matching pouch

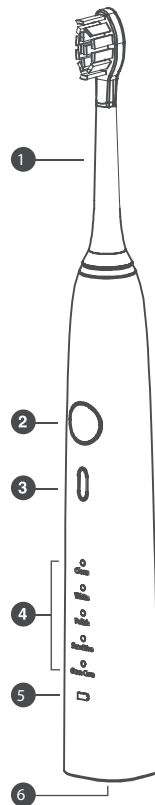The NEOSONIC® and the kids brush heads fit both the NEOKIDS® and the NEOSONIC® toothbrush handles. (Use sensitive or gum care mode for NEOSONIC® handle)

And don't forget to complete your brushing by **cleaning between your teeth with NEOPICK®, our dental floss with handle, or with one of our NEOPULSE® dental water jets**, because a tooth has not just 3 sides but 5!

Your NEOSONIC® toothbrush is safe for braces and restorations (crowns, dental veneers ...), we recommend that you use the «Ultra Soft» brush head. With orthodontic appliances, the brush heads may wear out more quickly. We advise you to change them every 2 months.

MAINTENANCE ADVICE

After each use, remove all toothpaste residues:

- Step 1: remove the brush head and rinse the bristles well with warm water
- Step 2: rinse the metal shaft area with lukewarm water
- Step 3: clean the handle with a damp cloth
- Step 4: **dry the whole area with a soft, clean towel** and put back the protective cap.

Maintenance and disposal:

If the toothbrush is to be stored for an extended period of time (six months or more), be sure to fully charge it before storing it in a cool, dry place. This device is equipped with a lithium-ion battery. When the battery reaches the end of its life, recycle the toothbrush at the nearest recycling centre in accordance with local requirements. The packaging box must be placed in the paper recycling bin. You can also contact us via the « assistance » form on our website (www.neopulse.fr) so that we can give your device a second life. In order to extend the life of our NEOSONIC® toothbrushes, NEOPULSE® has a repair service in (Yvelines, France).

2 YEAR WARRANTY

NEOPULSE® warrants that this product is free from material and manufacturing defects for a period of two years from the date of purchase. Keep your purchase receipt as proof of the date of purchase.

NEOPULSE® reserves the right to request that the product be returned for analysis to confirm that the device is covered by the warranty. NEOPULSE® will repair or replace the product, provided that the product has not been misused or abused, that it has not undergone any modification or damage after its purchase and that it has been used in accordance with these instructions. For any repair or replacement service offered under the warranty, please contact Neoteam Distribution by email at conso@neopulse.fr

Upon customer request, a copy of the warranty will be provided by email. This limited warranty does not cover accessories and consumable parts such as brush heads, USB-C cable, etc.

FIRST USE:

How to use the toothbrush?

1. Mandatory, charging the device before the first use


Plug the USB-C cable into a universal AC adapter and place the power cord into the charger plug (6) of the device, for 24 hours for the first use (even if your toothbrush arrives already charged). The battery level indicator (5) on the device changes from flashing green to a fixed green when the toothbrush is fully charged.

Do not operate the toothbrush while it is in charge. Don't forget to unplug it to use it. Under normal conditions of use (1 person brushing twice a day for 2 minutes), a full charge allows 3 months of use. When the battery needs to be recharged, the battery level indicator (5) will flash red. A full charge takes about 4 hours.

2. Use of the toothbrush:

Step 1: Brush head set up

Align the brush head so that the bristles are in line with the front of the handle, then slide it along the shaft and clip it in place.

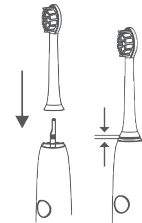To remove it: pull in the direction of the brush and, above all, do not shearing, as this could damage the brush motor (the warranty will be invalid).

Step 2: Customize your toothbrush by choosing from 15 possible combinations: 5 brushing modes X 3 strength settings

Press the setting button (2) and select the brushing mode (4), that means the **vibrations amplitude** that best suits your teeth and gums:

- **Clean:** standard cleaning mode to effectively remove plaque
- **White:** powerful cleaning mode for deep cleaning
- **Polish:** strong vibration mode that alternates between clean mode and white mode, ideal for removing surface stains and massaging the gums
- **Sensitive:** gentle brushing mode
- **Gum Care:** mode for very sensitive gums and teeth

Press the setting button (2) several times within 3 seconds to switch from one mode to another.


Then select the strength (3), that means the **power of the vibrations**:

- **Bright white:** high
- **Medium white:** medium
- **Low white:** low

With the toothbrush switched on, press button (3) to change the strength setting. Always start with the lowest setting and then find the strength that best suits your needs.

Once your toothbrush has been customized to one of the 15 settings, it is stored in memory. Just press the ON button (2) to start brushing.

neopulse

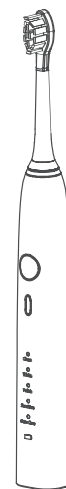www.neopulse.fr