Swop Physio

Instruction Manual

Treatment Tables for Physiotherapy

SWOP5.2.A66.B.PRO
SWOP5L.2.A66.B.PRO
SWOP2L.1.A66.B.EVO
SWOP4.1.A66.B.EVO
SWOP5.2.A66.B.EVO
SWOP7.2.A66.B.EVO

SCAN HERE
to read the
manual in
other
languages

VERSION 8.2026

Instruction Manual

Treatment Table for Physiotherapy

Table of contents

- ♥ 1. Introduction
-  2. Use of the manual
-  3. Warranty
-  4. Product Liability
-  5. EC Declaration of conformity
-  6. Safety instructions
-  7. Dimensions
-  8. Features & Operation
 - Wheel Block/Unblock system
 - Headrest adjustment
 - Breather hole and plug
 - Handlebar support
 - TwinLift
 - Peripheral bar
 - Switches for lumbar flexion
 - Backrest adjustment
 - Paper roll holder
 - Side widening supports
 - Safety lock
 - Interchangeable front panels
 - LED lighting
 - Bamboo sound system
 - Heating
-  9. Maintenance and cleaning
-  10. Troubleshooting
-  11. Service
-  12. EMC Guidance

1. Introduction

CONGRATULATIONS! You took NAGGURA home!

Thank you for having acquired a NAGGURA product and for the trust shown in us.

We are sure that this medical examination table will turn into the pillar of your work and will more than fulfil your most demanding work requirements. NAGGURA's experts are in the lead of a project created with great devotion. Every single detail has been cared to join technology, performance and design altogether in a single table. Our aim is designing and manufacturing INTUITIVE, DYNAMIC and HIGH-TECH tables of a DARING and UNIQUE design.

Thanks for your choice and start enjoying the NAGGURA experience.

 naggura'

2. Use of the Manual

Before using the treatment table, please read the instruction manual carefully. It is important to keep this manual and store it in a safe place so that it can be consulted at any time.

By carefully following these instructions you will be able to avoid any problems in using the treatment table.

Some texts and/or illustrations may not belong to your Naggura table since they are general and it depends on your product choice at the purchase time.

3. Warranty

The warranty period is effective on the date of purchase of the treatment table for a period of two (2) years. It covers all the components of the treatment table under the normal and correct use. It also covers defects in workmanship from factory. The above-mentioned warranty does not cover defects or malfunction caused by:

- Misuse, abuse or improper maintenance of the table and its upholstery;
- Non-observance of the precautions and warnings stated in this instruction manual;
- Lack of maintenance or repair tasks carried out by non-qualified service personnel who are unauthorized by the manufacturer;
- Wrong and faulty electrical installations or which do not fulfil the legal requirements;
- Naggura company is not liable for incidental or consequential damage to property or business.

4. Product Liability

Under any circumstances, the producer shall be liable for any indirect, special, incidental or consequential damage incurred by the user, the patient or the table itself as a result of:

- a. A wrong diagnosis;
- b. Incorrect use of the table or any of its accessories;
- c. The wrong interpretation or non-observance of this user manual;
- d. A lack of maintenance of the table;
- e. Or maintenance works of the table and/or repairs performed by persons unauthorized to do so by the manufacturer.

According to the current legislation, the manufacturer cannot be held responsible for any flaws in the table or its components and accessories after a period of ten (10) years from the date on which the producer puts into circulation the actual product. Therefore, producer cannot be obliged to supply any spare parts or accessories for this table.

Expected lifespan

Provided that the operator and user takes care of the table and undertakes regular maintenance tasks as described in section 9 of the current manual, the expected lifespan of the Naggura table can reach over 10 years from production date.

The intended use of the Treatment Table is to position a patient for the treatment. The unit is intended for use by qualified professional within a professional setting.

5. EC Declaration of conformity

Naggura products are produced and marked in conformity with the directives hereby stated :

Low voltage : 2014/35/CE

Electromagnetic compatibility : 2014/30/CE

These products are manufactured according to the essential requirements established in the European Medical Devices Directive (Annex II+III of Regulation (EU) 2017/745):

EN ISO 14971: 2019/A11:2021, EN ISO 15223-1:2021/A1:2025, EN ISO 20417:2021,
EN ISO 10993-1: 2020, EN ISO 10993-5: 2009/A11:2025, EN ISO 10993-10: 2023, EN ISO 10993-23: 2021, ISO 13485:2016,
EN 60601-2-52:2010+A1:2015, EN 60601-1: 2006/A13:2024, EN 60601-1-2: 2015+ A1:2021, EN 62366-1: 2015/A1:2021.

The Swop line has a guarantee of two (2) years in all parts. The manufacturer will not be held responsible for the results of maintenance or repairs performed by unauthorised persons. For further information on guarantee terms, please contact your distributor.

6. Safety instructions

Responsibility of the owner

Naggura treatment tables are only meant to be used by professional personnel in physiotherapy, rehabilitation, examination, massage and aesthetics fields.

It is the responsibility of the owner to ensure that training is given to all operators and third parties (like for example cleaning personnel) on safe use of the treatment table, and that they have read and understood all the operating and safety instructions. Only trained operators may use the treatment bed. Please keep the manual in a place where it is easily accessible to all users, and make it an important part of your employee training program.

Safety on installation and electricity


- Place the table in a space with free access from the 4 sides. Do not place the table against the wall or against any obstacles.
- Do not use the table in a high-humidity environment or in so-called "wet rooms" (hydrotherapy spaces). Any prolonged exposure to high moisture conditions may affect the electrical components and may result in failure.
- Check the technical data and electrical parameters to ensure that the supply voltage meets the technical specifications of the table. If not, contact your supplier before connecting the device.
- Make sure that the unit is correctly plugged into an earth intake, following the local and national requirements and standards.
- Make sure that the power cable is correctly and safely placed in such a way that it cannot get trapped between any mechanism.
- Pay special attention that the supply cable does not get pinched under the feet of the table, (e.g. after having moved the table on its wheels).
- If the table accidentally gets wet, unplug it immediately.
- Do not leave children unattended near the table.
- Unless being monitored or instructed by personnel who can be responsible for their personal safety, the people who have physical obstacles, tactile disorders, mental disorders, lack of experience, common sense or children, are not allowed to use the equipment.
- Do not use the device near electromagnetic transmitters. Read and follow all the information about electromagnetic compatibility (EMC) in page 21.
- If the power cord or plug is damaged, please contact the distributor, service center or experienced technicians for repairing.
- Please pay attention that the continuous working times should be "Max.2Min. ON; 18Min. OFF", when operating the motor.

Safety instructions

Safety during treatments

- Don't let anybody, except the operator and the patient, in the treatment room.
- The table is not intended to transport patients. Therefore it cannot be used for this purpose.
- Do not sit on the head section or end of the leg section of the table. The head section and armrests are not suitable to sit on.
- Patients are not to be left unattended,
- Never leave the table on its wheels during the treatment.
- Do not exceed the maximum load capacity of the table.
- Place the table on a stable and horizontal surface. Level it if necessary with the levelling feet under the table.
- Guide the backrest section with the other hand when adjusting the backrest by gas levers.
- Ensure that neither the patient nor the operator can get any part of the body or clothes trapped between the moving parts of the mechanism when the table is being raised or lowered.

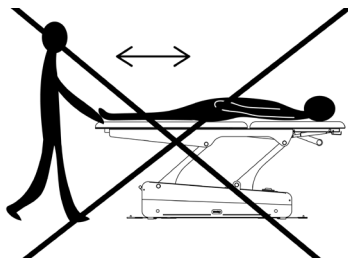
Safety before or after the treatment

- Do not leave anybody unattended alone near the table.
- Put the safety lock in "OFF" position and remove the key.
- Disconnect the table if not used for a long period of time.

Safety when cleaning the bed or the therapy room


- Disconnect the table from the power supply.
- Do not get under the table or between its moving parts.

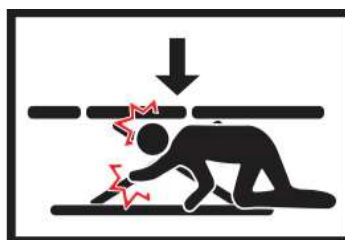
7. Product dimensions

Swop 5L Pro

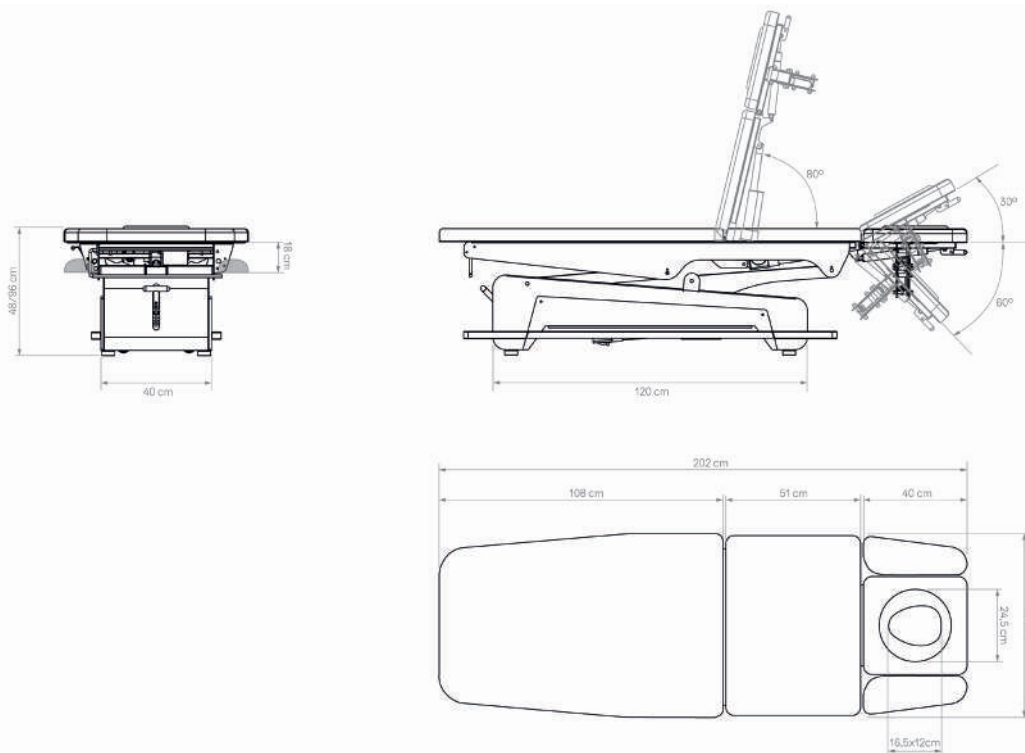
Product size	202 (+/- 37,5 cm with movement) x 70 x 48,5 / 96 cm
Package size	194 x 74 x 54 cm
Max loading weight	200 Kg
Net weight	109 Kg
Gross weight	129 Kg
UDI	0 697459027 095 9
Heating	NO
Upholstery type	Naggura Skin PU in Grey
Upholstery density	wooden base (1,5cm) + foam 40 kg/m ³ (3,5cm) Total thickness: 5cm
Accessories	A66.COVER.SWOP5, MUSIC.SWOP, SW0049

Swop 5 Pro

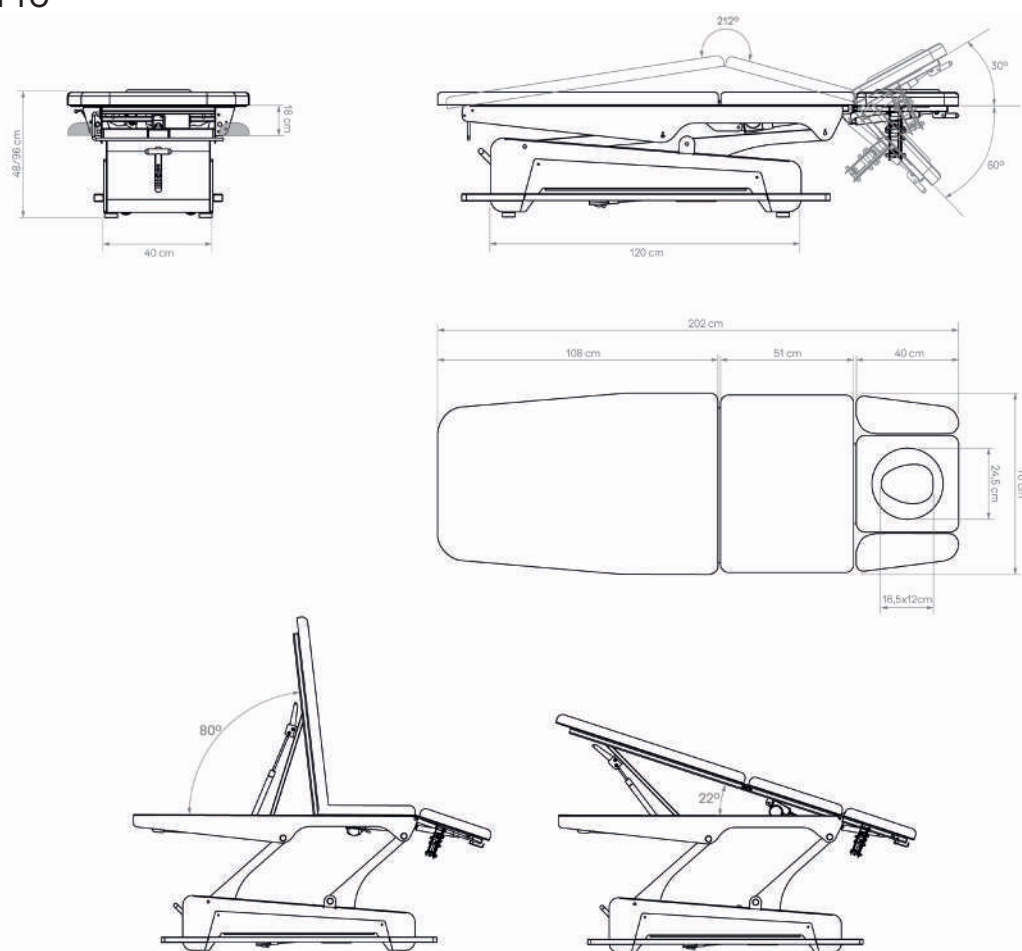
Product size	202 (+/- 37,5 cm with movement) x 70 x 48,5 / 96 cm
Package size	194 x 74 x 54 cm
Max loading weight	200 Kg
Net weight	109 Kg
Gross weight	129 Kg
UDI	0 697459027 094 2
Heating	NO
Upholstery type	Naggura Skin PU in Grey
Upholstery density	wooden base (1,5cm) + foam 40 kg/m ³ (3,5cm) Total thickness: 5cm
Accessories	A66.COVER.SWOP5, MUSIC.SWOP, SW0049

Swop 2L EVO

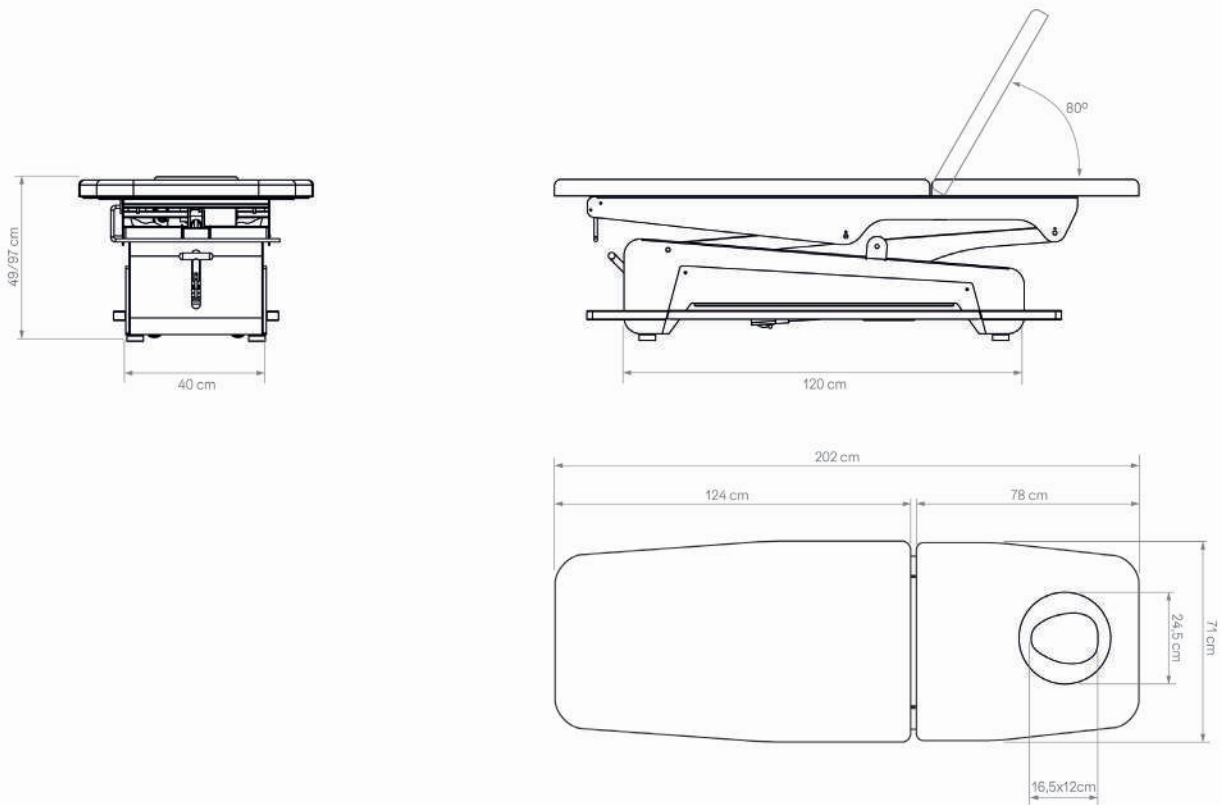
Product size	202 (+/- 37,5 cm with movement) x 71 x 49,5 / 97 cm
Package size	206 x 74 x 55 cm
Max loading weight	200 Kg
Net weight	100 Kg
Gross weight	120 Kg
UDI	0 697459027 096 6
Heating	YES
Upholstery type	Naggura Skin PU in Grey
Upholstery density	wooden base (1,5cm) + foam 60 kg/m ³ (4cm) + 40 kg/m ³ (1.5cm) Total thickness: 7cm
Accesories	COVER.PHYSIO, MUSIC.SWOP

Swop 4 EVO

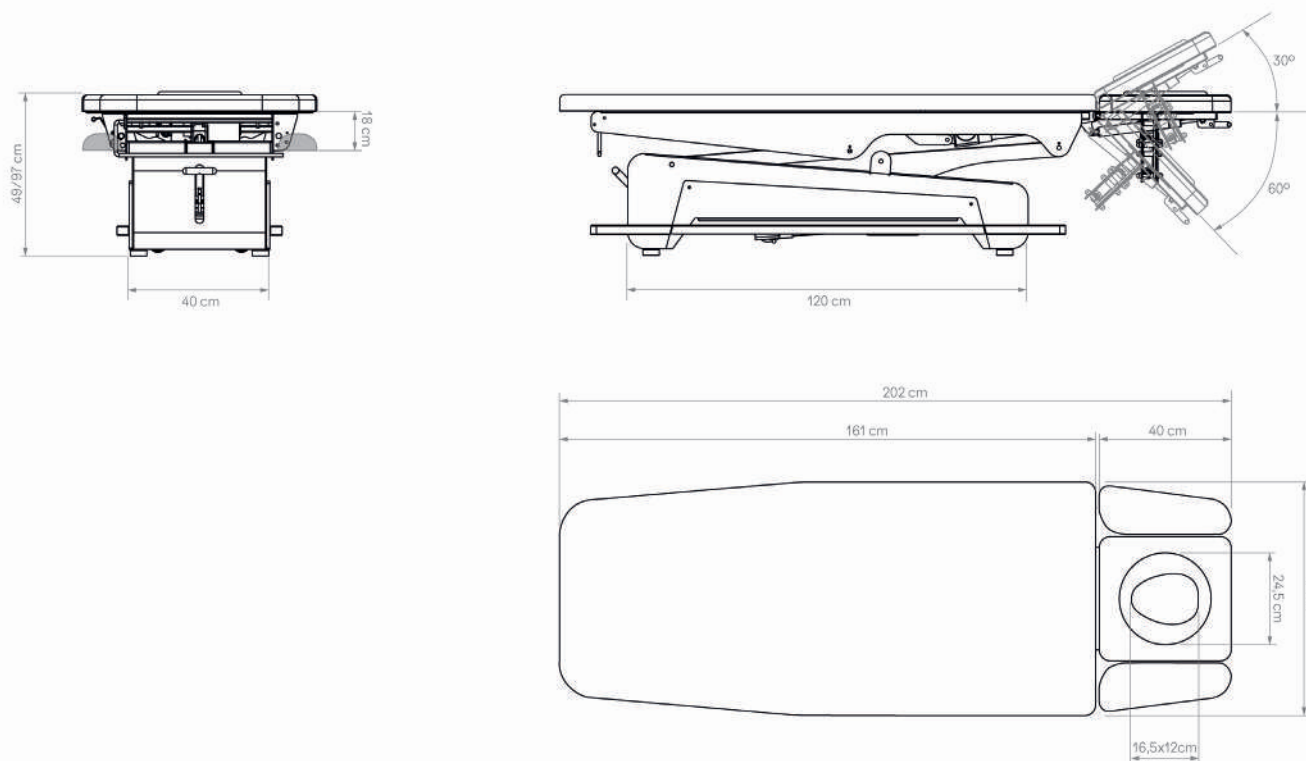
Product size	202 (+/- 37,5 cm with movement) x 71 x 49,5 / 97 cm
Package size	194 x 74 x 55 cm
Max loading weight	200 Kg
Net weight	100 Kg
Gross weight	120 Kg
UDI	0 697459027 097 3
Heating	YES
Upholstery type	Naggura Skin PU in Grey
Upholstery density	wooden base (1,5cm) + foam 60 kg/m ³ (4cm) + 40 kg/m ³ (1.5cm) Total thickness: 7cm
Accessories	A66.COVER.SWOP5, MUSIC.SWOP

Swop 5 EVO

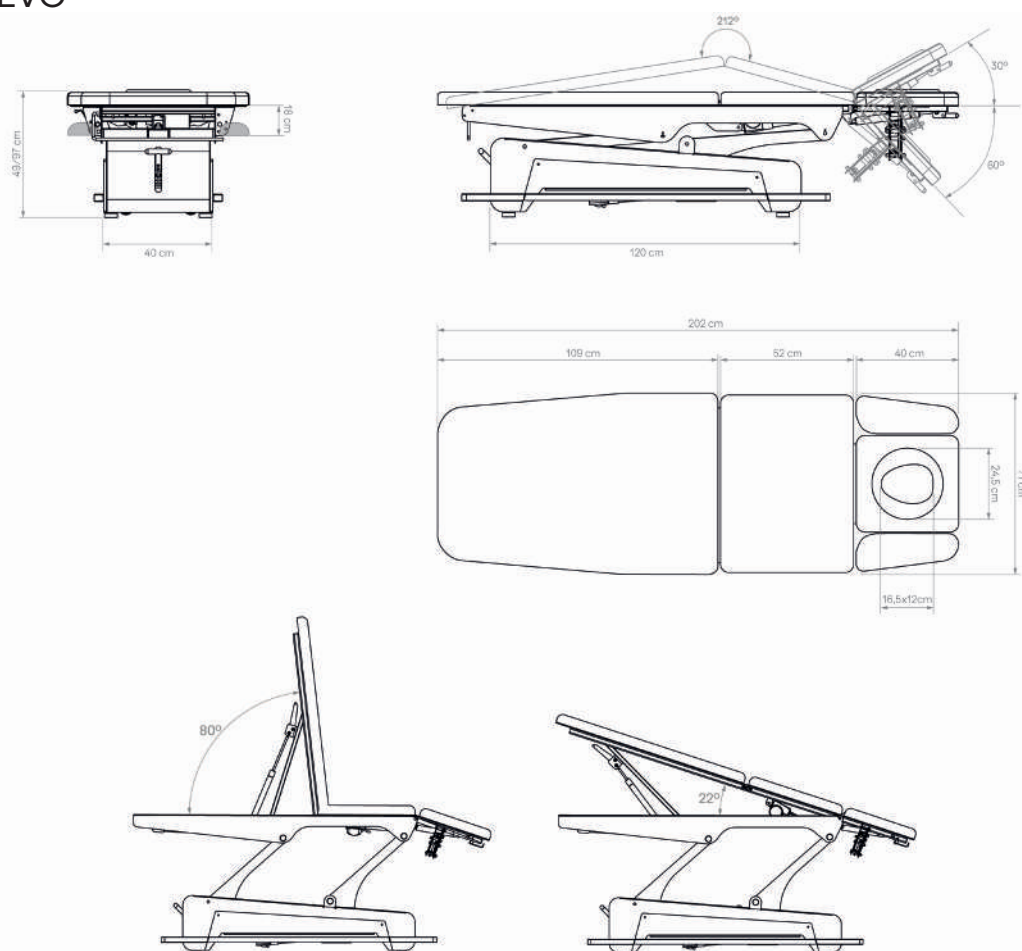Product size	202 (+/- 37,5 cm with movement) x 71 x 49,5 / 97 cm
Package size	194 x 74 x 55 cm
Max loading weight	200 Kg
Net weight	125 Kg
Gross weight	145 Kg
UDI	0 697459027 098 0
Heating	YES
Upholstery type	Naggura Skin PU in Grey
Upholstery density	wooden base (1,5cm) + foam 60 kg/m ³ (4cm) + 40 kg/m ³ (1.5cm) Total thickness: 7cm
Accessories	A66.COVER.SWOP5, MUSIC.SWOP, SW0050


Swop 5S EVO

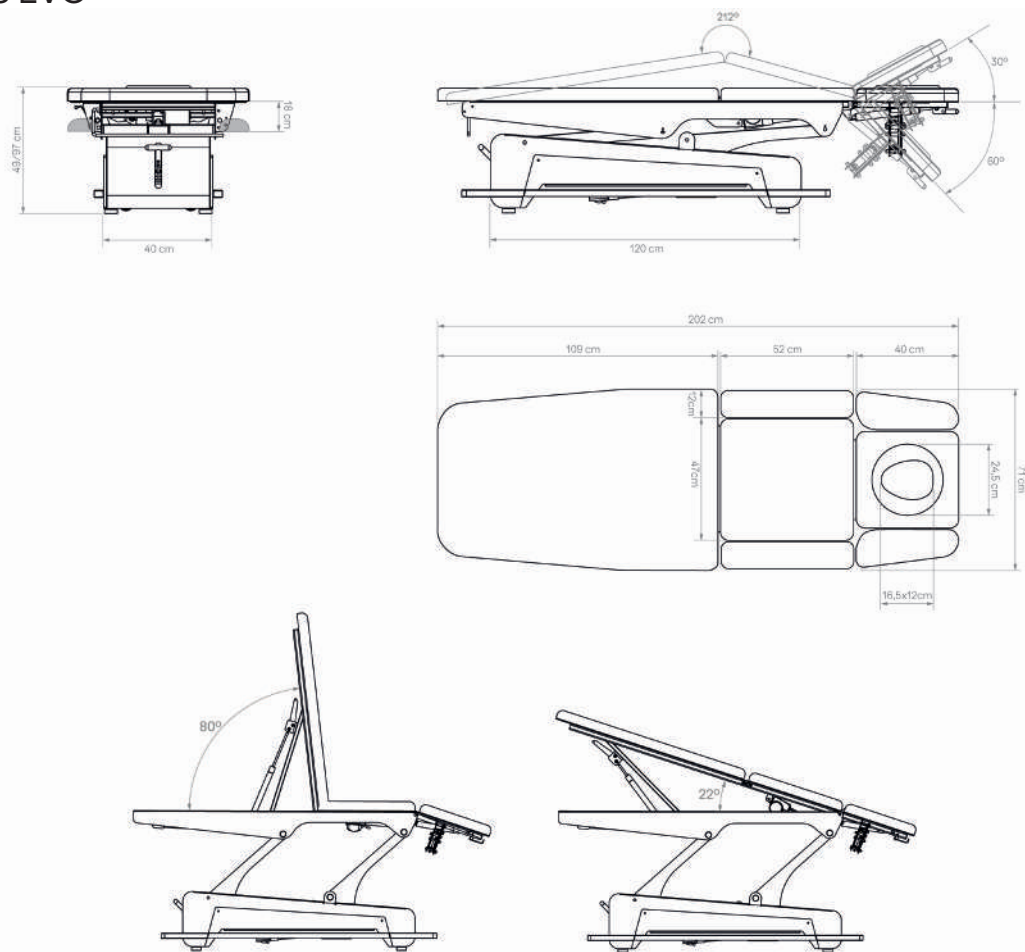Product size	202 (+/- 37,5 cm with movement) x 71 x 49,5 / 97 cm
Package size	194 x 74 x 55 cm
Max loading weight	200 Kg
Net weight	125 Kg
Gross weight	145 Kg
UDI	0 697459027 099 7
Heating	YES
Upholstery type	Naggura Skin PU in Grey
Upholstery density	wooden base (1,5cm) + foam 60 kg/m ³ (4cm) + 40 kg/m ³ (1.5cm) Total thickness: 7cm
Accessories	A66.COVER.SWOP5, MUSIC.SWOP

Handling and transport

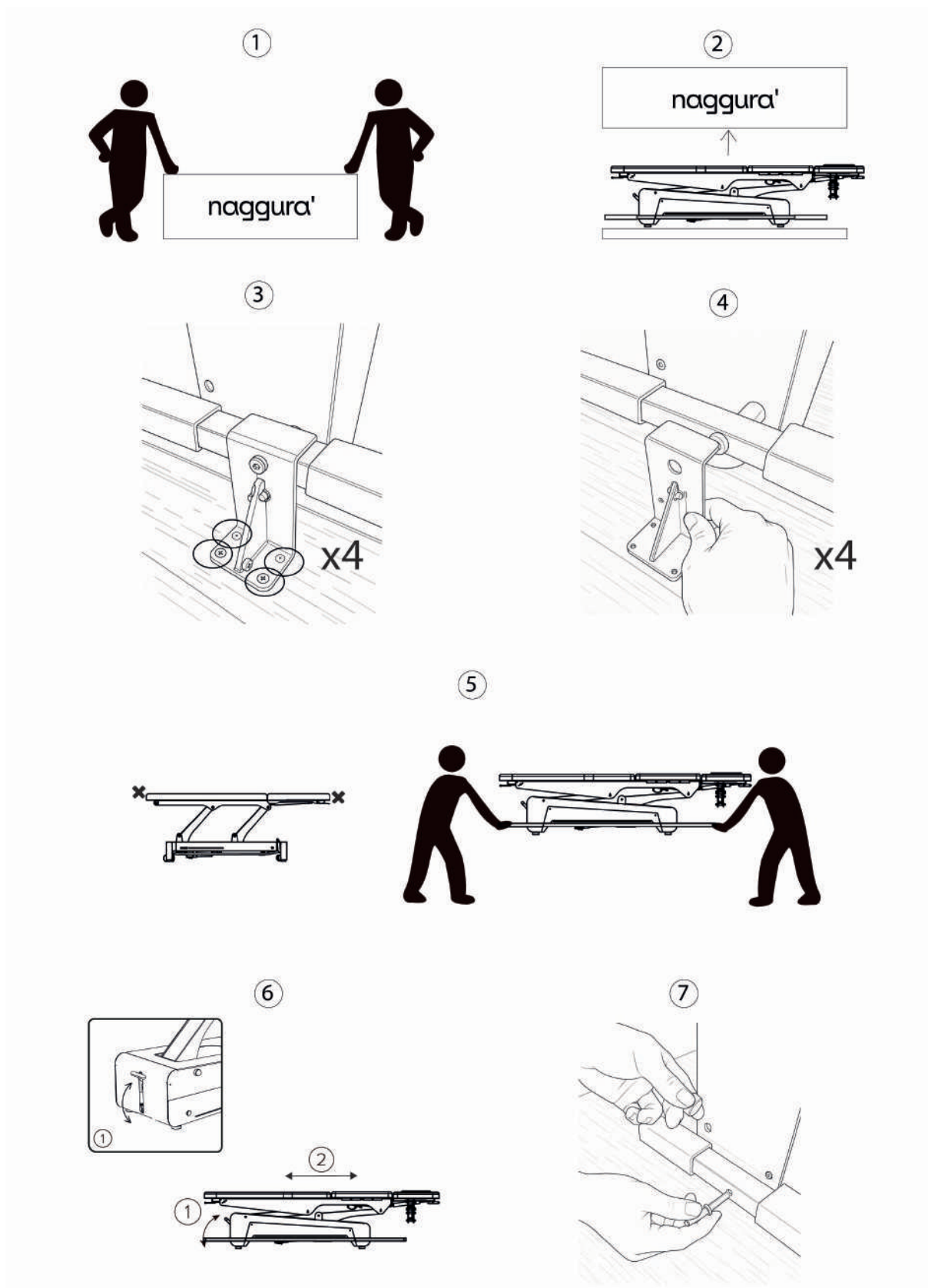Specifications


Power	AC 110 - 240V, 50/60Hz
Max power	300W
Mains cable length	250 cm
Environment temperature	0°C - +40°C
Relative humidity	30% - 90%
Loading capacity	160 Kg
Electrical class	I
Applied parts	Type B
Medical Classification	Class I

8. Operation

Power input


Voltage and frequency input shall be 220-240V, 50/60Hz. Plug the table cable into the socket as shown below:





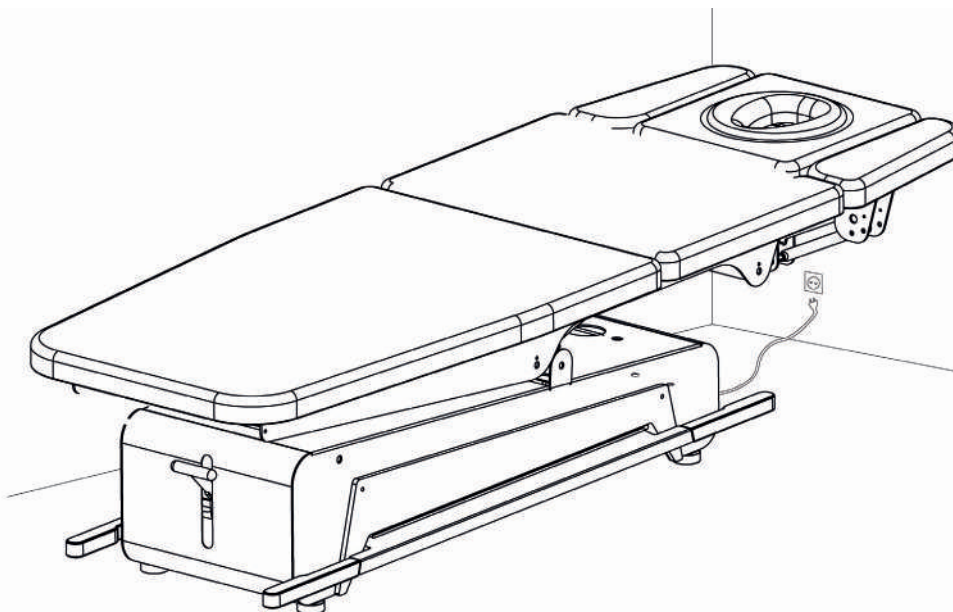
Operation of Swop Physio

1. Wheel block/unblock system
 2. Lever to adjust the headrest
 3. Breather hole and plug
 4. Paper roll Holder
 5. TwinLift
 6. Peripheral switch
 7. Switches for upper motor
 8. Switches for heating
 9. Lever to adjust the backrest
 10. Side widening supports
 11. Safety lock
 12. Interchangeable front panels
 13. LED lighting
 14. Bamboo sound system
- Optional

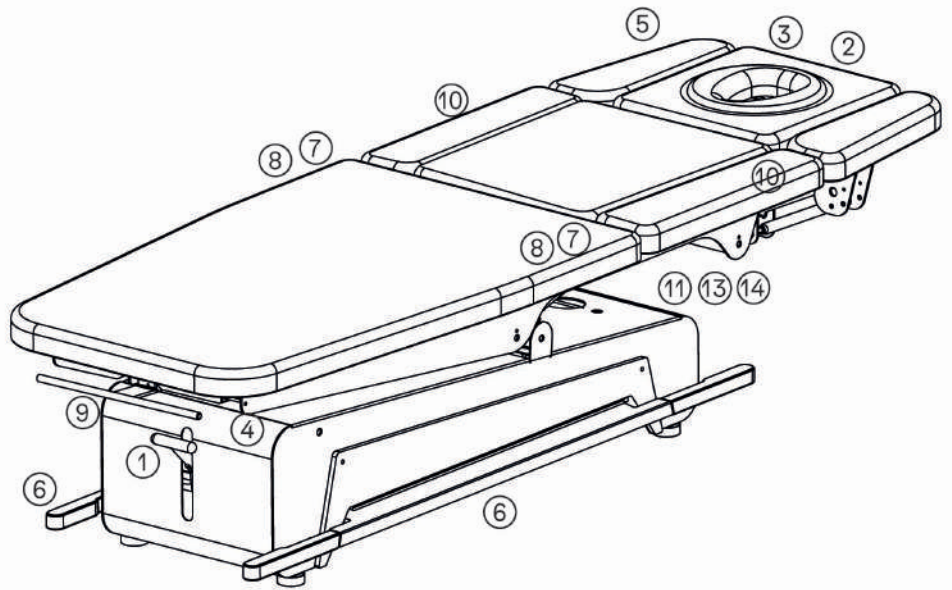*Drawings might not be the same as your table

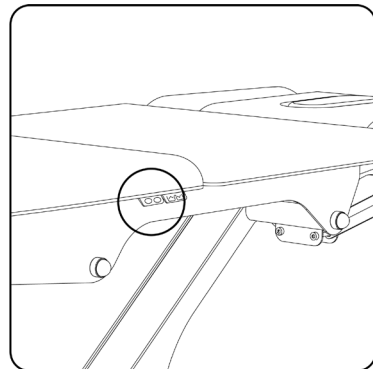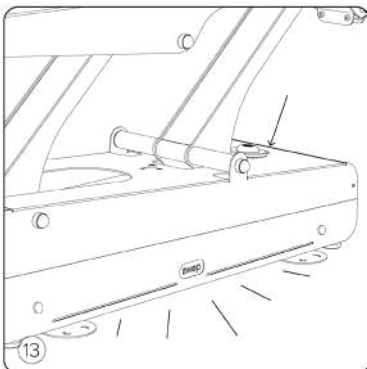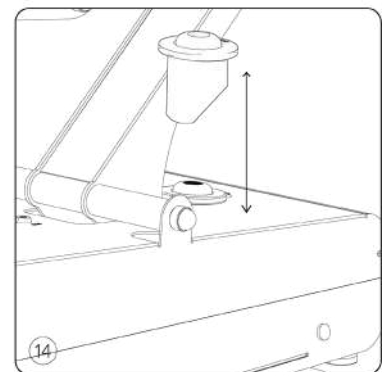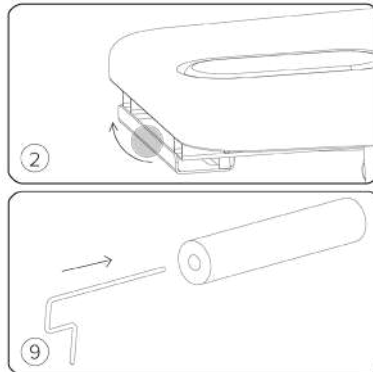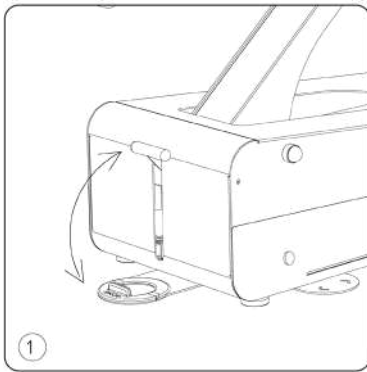
Operation


Peripheral Bar - How to adjust height

The height of the table can be adjusted using the peripheral bar equipped as standard.

To lift the table, simply press down the side bar with your foot. Once the desired height is achieved, release the button.
To lower the table, push up the side bar, until the desired height.

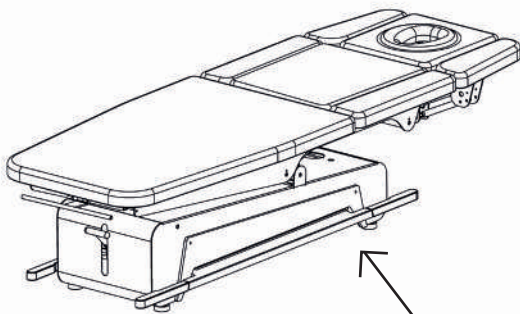
PRESS DOWN = RISE HEIGHT

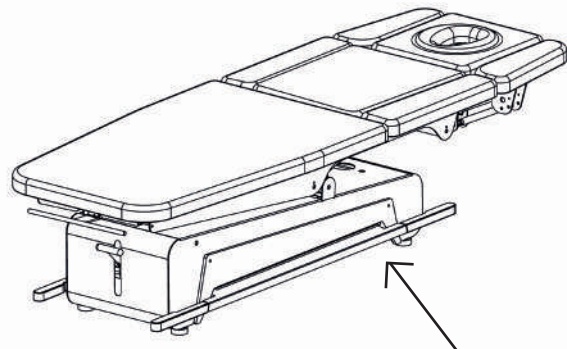PUSH UP = LOWER HEIGHT

The length of the peripheral bar can be regulated.
Simply extend each of the extremes pushing with your hands to the desired length.

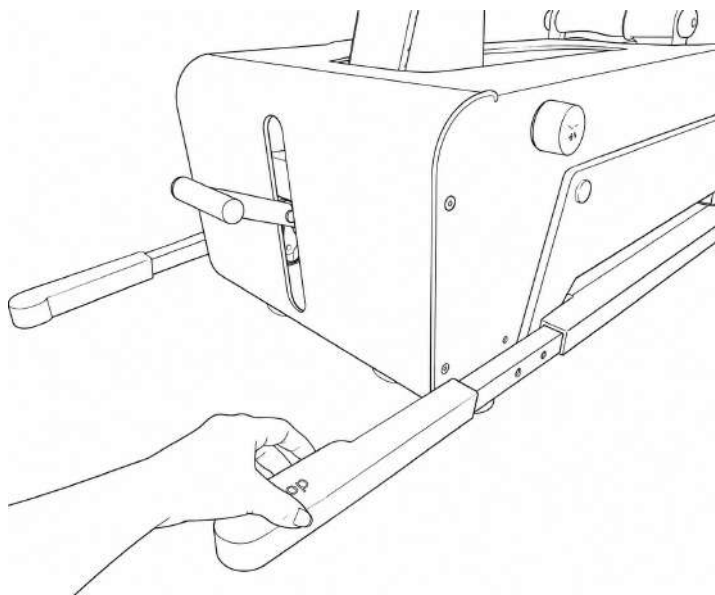
Operation


Switches for tables with 2 motors

If your Swop is a 2-engine table (Swop 5 Pro, Swop 5L Pro, Swop 5 EVO or Swop 5S EVO), the middle section (lumbar flexion) is adjusted by the push button box placed at both sides of the table.

In Swop 5L PRO, instead of lumbar flexion, the middle section + headrest are raised in chair position.

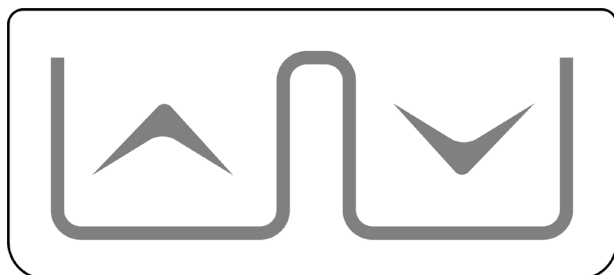Swop 5 Pro, Swop 5L Pro, Swop 5 EVO and Swop 5S EVO button panel for controlling second engine.

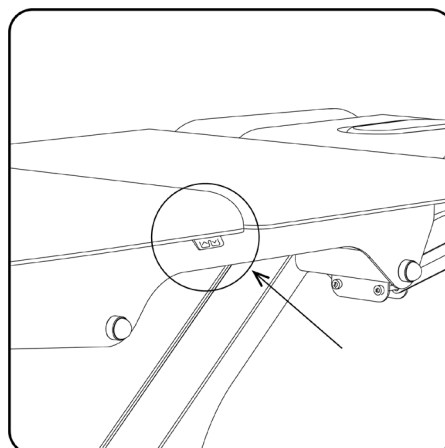
Switches for heating

One of the Swop EVO main characteristics is the integrated two-level heating function to warm up the cushions of the table.


It can be activated and regulated using the switches found on both sides of the table.

Press the symbol on the left to activate low intensity heating and press the symbol on the right to activate high intensity heating. To turn it off, press the activated button again.

The heating will automatically turn off after 1 h from the moment of activation.

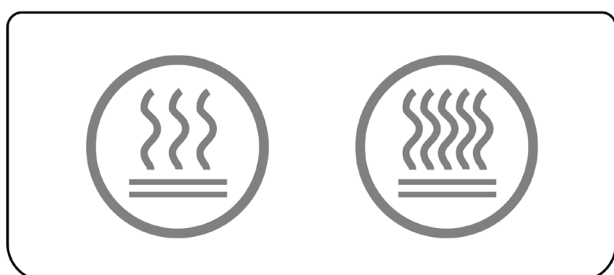Swop EVO button panel for controlling the heating function.

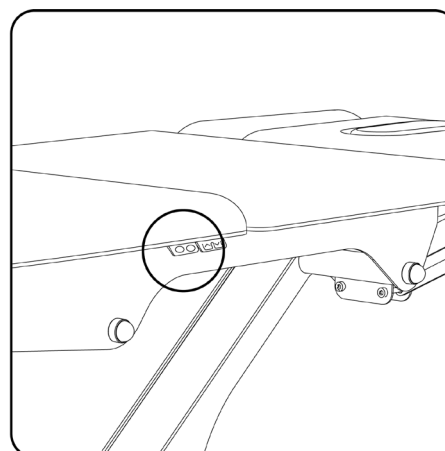
Operation

TwinLift Armrests

Swop 5 Pro, Swop 5L Pro, Swop 4 EVO, Swop 5 EVO and Swop 5S EVO are equipped with a head section with built-in arm supports. You can adjust them as it is shown in the pictures below. Press the lever on the side and bring the arm rests up and down to the desired position. They will move at the same time effortlessly. The patients can adjust them by themselves.

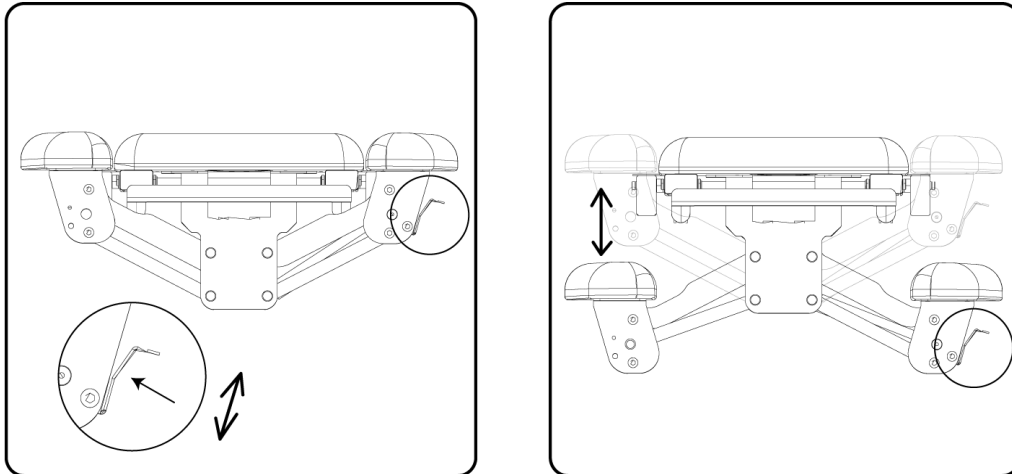
Safety lock

The table can be locked or unlocked using the key on the base. If the key is switched to ON-position, the table and all its functions are unlocked and ready to use.

If the key is in OFF-position (horizontal position), the table is locked and all the electrical systems are disabled.

The safety-key can only be removed when the safety lock is in OFF-position.

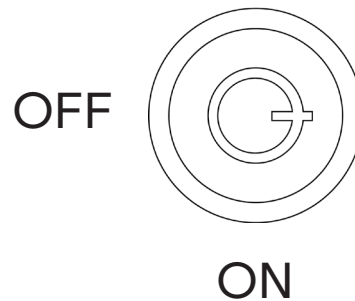Integrated safety system

All SWOP physio models have 2 integrated safety systems:

1-Double tip safety system: Electronic configuration mandatory in some EU countries.

Press 10 consecutive times the peripheral bar to activate the “double tip safety system”.

After a “beep” sound, the system will be activated you'll need to press two times the peripheral bar in order to raise or lower the height of the bed.

To deactivate, press again 10 consecutive times the peripheral bar.

2-Clutch motor: height motor with an automatic clutch system. In case somebody or something gets stuck between the upper and lower frame of the bed (when the bed moves downwards), the motor clutch will make sure that the motor stops pulling downwards. Only the “passive weight” of the bed will impact the object and there will be no additional pulling force of the motor.

Important: in case somebody or something gets stuck, make sure that 2 persons lift up the upper frame of the bed to release the object. After the object is fully released and safe, gently release the upper frame until the bed is in its lowest position.

Operation


Interchangeable Front Panels

The Swop treatment tables highlights by a very special feature: the front exchangeable panels. There are lot of different possibilities to dress your table up. Get a new table in a matter of seconds and according to your taste, the interior design of your centre, your mood or however you feel today!

The exchangeable panels are 2 units: one for one side of the table base and the other for its opposite.

The available panels are made of different materials (metal, wood and pmma) and multiple colours. You can visit our website to discover the full range of possibilities.

Changing the panels is so easy! Simply unscrew the 2 screws that subject the panel with your fingers (no tools required). Remove the current panel to place a different one.

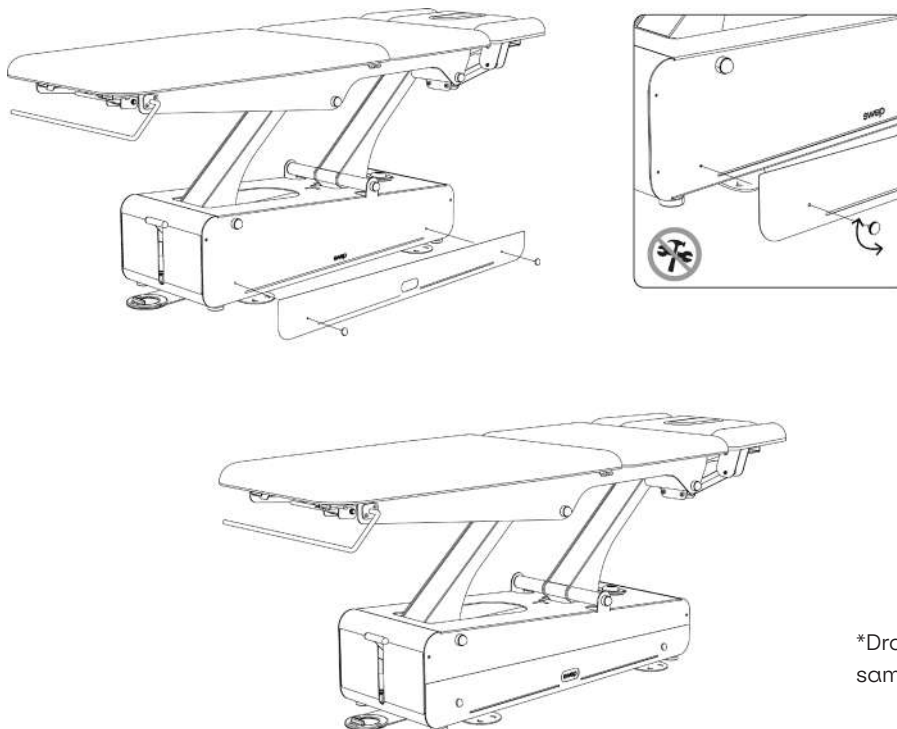*Drawings might not be the same as your table


Bamboo sound system

If you purchased the Bamboo sound system, your swop beauty bed arrives with an awesome Bluetooth speaker: we call it the Bamboo sound system.

The speaker is lodged inside a bamboo case with bevelled base. You can place the speaker on the table housing designed for this purpose or take it off and place it on your working table. Thanks to its inclined base you can direct the sound wherever you want to.

The Bluetooth speaker has a USB charger and additional instructions. Please read them carefully.

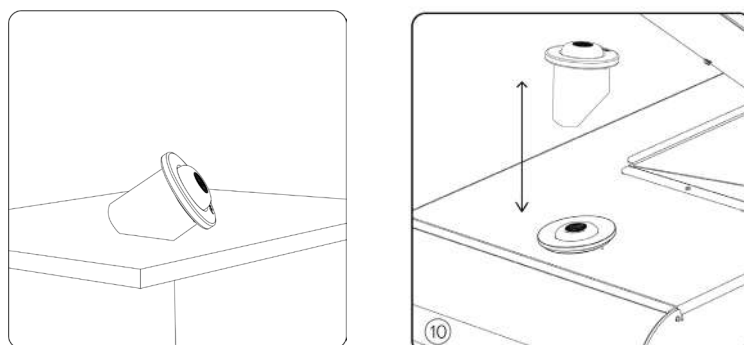9. Maintenance and cleaning

Cleaning the chassis and metal frame:

Clean the table with a clean, soft and dry cloth.

For persisting stains, use a white cloth, damp it with clear water and a domestic household cleaner FREE of alcohol.

Cleaning the cushions:

The Naggura tables have a high quality vinyl upholstery protecting against stains and abrasion.

The vinyl should be cleaned periodically in order to maintain its appearance and prevent build-up of dirt and contaminants. Any stains should be removed immediately.

- Cleaning after each patient:

If disinfection after each patient is required, use one of the products "list of disinfectants". In case of spilling massage oil, remove it with a PH-neutral soap immediately (1:9 liquid soap, neutral PH and water). Always remove residu with clean water after cleaning with soap.

- Cleaning at the end of each day:

In order to prevent the build-up of contaminants and grease, it is important that the surface is cleaned every day. Clean the entire surface using a mixture of 1:9 liquid soap (neutral PH) and water. Do not use soap or cleaning solutions that contain alcohol, ketones, xylene, acetates or solvents. Always remove residu with clean water after cleaning with soap.

- Cleaning each month or in case of stains:

Use the original Naggura Cuirnet high-performance cleaner in a dilution of 1:9 Cuirnet and water. If the stain cannot be removed, the cuirnet cleaner maybe be used undiluted.

- After the cleaning, always rinse with clear water and dry.

- NEVER use the Naggura Cuirnet cleaner for the metal frame of the table.

Maintenance

Naggura tables have been designed for an intensive use during long time thanks to its stability, reliability and durability. However, we advise to keep the table in perfect conditions and make a correct use of it. Unplug the table if you will not use it for a long period of time.

Environment information

If this symbol appears on electrical products or on its packaging, this means that they must be recycled properly. Please ensure that you are well informed of the local rules and regulations regarding to the removal of equipment and accessories.

10. Troubleshooting and service

The backrest moves down alone

The supporting gas lifts are loose.

- In such a case, remove the gas lift from the lower part and loosen it (unscrew it) 2 turns and fix it again. Proceed the same way with the second gas lift.

The back rest is blocked (it can't be lowered nor raised up)


The gas lift activating pin does not activate the lever. Remove the lower side of the gas lift and tighten it (screw it 2 turns until you see that the pin can activate the lever. Proceed the same operation with the second gas lift.

The up and down motor does not work

- Check if the power cable is disconnected, either from the wall socket or the control box.
 - Check the power cable for any damage.
 - Check all the plugs and connections.
 - There are 3 possible components that may be faulty: the peripheral switch, the control box or the motor itself.
- Contact your dealer to proceed accordingly.

Squeaking noise when moving the table up & down (mechanical noise)


- The axes are adjusted and greased from factory but with the time and use, they may need some little adjustment. Remove the grey plugs of the 2 connecting rods (4 at each side, so 8 in total). Loosen the axes a little bit by unscrewing them (approx. 2 turns). Apply a bit of neutral grease and plug them again.

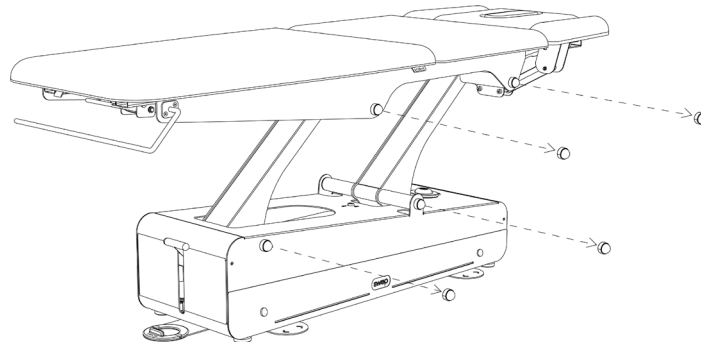Persisting problems

If the problem persists, contact your local distributor and provide the serial number and the model of the table, or a picture of the typeshield.

This can be found on the typeshield, located at the backside of the table, inside the top frame.

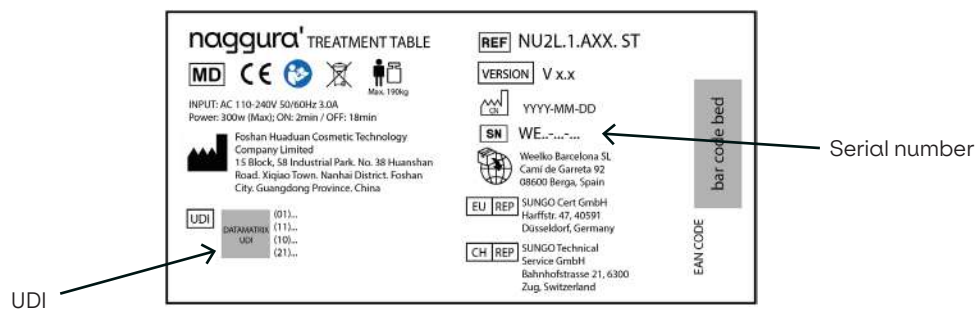11. Service

Only a technician authorised by Naggura is allowed to perform any kind of repairs on the bed.

If your technical staff is allowed by Naggura to repair the bed, Naggura can provide the necessary instructions and spare parts lists.

Please contact your local Naggura distributor for more information.

Manufacturer:

Foshan Huaduan Cosmetic Technology Company Limited
15 Block, 58 Industrial Park. No. 38 Huanshan Road. Xiqiao Town. Nanhai District. Foshan City. Guangdong
Province. China

Importer:

Weelko Barcelona SL
Camí de Garreta 92, 08600 Berga, Spain

EU REP:

SUNGO Cert GmbH
Harffstr. 47, 40591 Düsseldorf, Germany

CH REP:

SUNGO Technical Service GmbH
Bahnhofstrasse 21, 6300 Zug, Switzerland

12. EMC Guidance

Guidance and Manufacturer's Declaration - Electromagnetic Emissions	
Emissions test	Compliance
RF emissions CISPR 11	Group 1
RF emissions CISPR 11	Class A
Harmonic emissions IEC 61000-3-2	Class A
Voltage fluctuations/flicker emissions IEC 61000-3-3	Complies

Warning: It is not safe to operate the device in the presence of HF surgical equipment or places where the intensity of EM disturbances is high.

Warning: Use of this equipment adjacent to or stacked with other equipment should be avoided because it could result in improper operation. If such use is necessary, this device and other equipment should be observed to verify that they are operating normally.

Warning: Use only accessories and cables provided by the manufacturer. Accessories and cables not recommended may affect the EM emissions or immunity.

Warning: Maintain minimum distance of 12 inches (30 cm) between any part of the working table and portable RF communication equipment (including peripherals such as antenna cables and external antennas.) Performance of the table might degrade if proper distance is not maintained.

Guidance and Manufacturer's Declaration - Electromagnetic Immunity		
Immunity Test	IEC 60601-1-2 Test level	Compliance level
Electrostatic discharge (ESD) IEC 61000-4-2	±8 kV contact ±2 kV, ±4 kV, ±kV, ±15 kV air	±8 kV contact ±2 kV, ±4 kV, ±kV, ±15 kV air
Electrical fast transient /burst IEC 61000-4-4	±2 kV for power supply lines ±1 kV for input/output lines	±2 kV for power supply lines ±1 kV for input/output lines
Surge IEC61000-4-5	±1 kV line(s) to line(s) ±2 kV line(s) to earth	±1 kV line(s) to line(s) ±2 kV line(s) to earth
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	0% U; 0.5 cycle T at 0°, 45°, 90°, 135°, 180°, 225°, 270°, 315° 0% UT; 1 cycle 70% UT; 25/30 cycle 0% U; 250/300 cycle	0% U; 0.5 cycle T at 0°, 45°, 90°, 135°, 180°, 225°, 270°, 315° 0% UT; 1 cycle 70% UT; 25/30 cycle 0% U; 250/300 cycle
Power frequency (50/60 Hz) Magnetic field IEC61000-4-8	30 A/m 50Hz/60Hz	30 A/m 50Hz/60Hz
Conducted RF IEC61000-4-6	3 V r.m.s. 150 kHz to 80 MHz 6 V RMS in the ISM bands between 0.15 MHz and 80 MHz	3 V r.m.s. 150 kHz to 80 MHz 6 V RMS in the ISM bands between 0.15 MHz and 80 MHz
Radiated RF IEC61000-4-3	3 V/m 80 MHz - 2,7 GHz 80% AM at 1 kHz	3 V/m 80 MHz - 2,7 GHz 80% AM at 1 kHz
NOTE: UT is the a.c. mains voltage prior to application of the test level.		

Test Specifications for enclosure port immunity to RF wireless communications equipment

Immunity Test	IEC60601 test level				Compliance level
	Test frequency	Modulation	Maximum power	Immunity level	
Radiated RF IEC 61000-4-3	385 MHz	**Pulse Modulation: 18Hz	1.8W	27 V/m	27 V/m
	450 MHz	*FM+ 5Hz derivation: 1kHz sine	2W	28 V/m	28 V/m
	710 MHz 745 MHz 780 MHz	**Pulse Modulation: 217Hz	0.2W	9 V/m	9 V/m
	810 MHz 870 MHz 930 MHz	**Pulse Modulation: 18Hz	2W	28 V/m	28 V/m
	1720 MHz 1845 MHz 1970 MHz	**Pulse Modulation: 217Hz	2W	28 V/m	28 V/m
	2450 MHz	**Pulse Modulation: 217Hz	2W	28 V/m	28 V/m
	5240 MHz 5500 MHz 5785 MHz	**Pulse Modulation: 217Hz	0.2W	9 V/m	9 V/m
<p>Note* - As an alternative to FM modulation, 50% pulse modulation at 18 Hz may be used even if it does not represent actual modulation.</p> <p>Note** - The carrier shall be modulated using a 50% duty cycle square wave signal.</p>					

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